

Impressions



PrepQuick

Chromaclone

Chromaclone PVS

Clone Bite

Chromaclone PVS Accessories

Chromaclone PVS Mixing Tips



Dr. Robert Feasel, Shoreline, WA
Washington

PrepQuick™


PRE-IMPRESSION SURFACE CLEANING AGENT



Black Micro FX tip
#1434 - 500pk - pg 128

PrepQuick is a 2% glycolic acid, which removes contaminants from the preparation and minimally affects the smear layer. A clean preparation reduces surface tension and allows for a lower wetting angle, improving reproduction of detail and minimizing air bubbles. It does not interfere with the bond of temporary or permanent cements, and it rinses cleanly.

"I get beautiful impressions every time."
DR. LESLIE NULTY – SARASOTA, FL

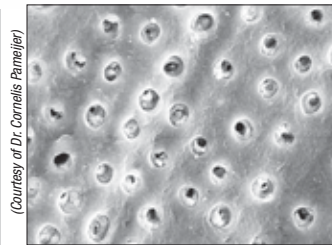

(Courtesy of Dr. Cornelis Farneljer)

PrepQuick cleans while leaving tubules occluded.



PrepQuick, rinse, dry, and impress.
Result: A quality impression.

Disposable Impression Trays



- Moldable
- Rigid to prevent distortion
- Provides excellent retention

#516 - PrepQuick Kit

4 - 1.2ml syringes
20 - Black Micro FX tips



Shelf Life: 24 months.

#517 - PrepQuick Refill

4 - 1.2ml syringes



Shelf Life: 24 months.

- | | |
|-------|---|
| #5070 | #5071 - Maxillary Left - Mandibular Right |
| #5072 | #5073 - Maxillary Right - Mandibular Left |
| #5074 | #5075 - Maxillary Large Size 1 |
| #5076 | #5077 - Maxillary Medium Size 3 |
| #5078 | #5079 - Maxillary Small Size 5 |
| #5080 | #5081 - Mandibular Large Size 2 |
| #5082 | #5083 - Mandibular Medium Size 4 |
| #5084 | #5085 - Mandibular Small Size 6 |
| #5086 | #5087 - Anterior |

12 - Trays

36 - Trays



Chromaclone®

COLOR-CHANGING ALGINATE

- Great Walterberry® and Bubble Gum flavors
- Color-changing
- Dust-free mixing
- Silicone-like consistency
- Unsurpassed accuracy

Chromaclone is a fast-setting, color-changing, Class A, Type I sodium alginate. It is formulated to be stiffer than conventional alginates to increase accuracy. The change in color indicates time for mixing, loading, and tray insertion. The quick setting time makes it perfect for gaggers.

"I've been a Chromaclone user for many years. We like the accuracy, even mixes, and quality results we get. The patients appreciate the flavor."
 DR. DAVID N. PETERSON – SAUK RAPIDS, MN



#389 - Water Cup and Scoop Set

1 - of each cup and scoop

Free of charge



Time in mouth

Regular Set - 80 seconds

Fast Set - 53 seconds

Super Fast Set - 45 seconds

#4032 - Walterberry Regular Set

1 - 1lb bag of Chromaclone
 1 - Plastic container, water cup, and scoop



Shelf Life: 60 months.

#671 - Bubble Gum Fast Set

1 - 1lb bag of Chromaclone
 1 - Plastic container, water cup, and scoop



Shelf Life: 60 months.

#5500 - Ortho Bubble Gum Super Fast Set

1 - 1lb bag of Chromaclone
 1 - Plastic container, water cup, and scoop



Shelf Life: 60 months.

#4031 - Walterberry Regular Set

1 - 1lb bag of Chromaclone



Shelf Life: 60 months.

#672 - Bubble Gum Fast Set

1 - 1lb bag of Chromaclone



Shelf Life: 60 months.

#5501 - Ortho Bubble Gum Super Fast Set

1 - 1lb bag of Chromaclone



Shelf Life: 60 months.

Chromaclone® PVS
HYDROPHILIC POLYVINYL SILOXANE
IMPRESSION MATERIAL



- Hydrophilic and thixotropic
- Dimensionally stable with great elastic recovery
- Allow you to pour at convenience
- Tear-resistant
- Accurate with multiple pours
- Time in mouth: Fast 2 minutes and Regular 3 minutes

Chromaclone PVS impression material is super hydrophilic not only during application, but also during immediate pouring of gypsum. It contains special reinforcing branched silicone resins for fabulous tear resistance, so impressions can be easily removed. The Bubble Gum flavor gives the patient a more pleasant impression-making experience. Chromaclone PVS remains dimensionally stable for weeks, allowing you to pour at your convenience.

The accuracy of any impression material, including Chromaclone PVS, can be increased with the use of Ultradent's Tissue Management System, pages 48-49. Profound hemostasis is paramount before making impressions.

#5000 - Chromaclone PVS Regular Set Kit
#5001 - Chromaclone PVS Fast Set Kit

- 1 - 50ml cartridge of each: light, medium, and heavy
- 1 - 10ml bottle tray adhesive
- 1 - Dispensing gun
- 19 - Mixing tips - variety
- 20 - IntraOral tips



Shelf Life: 36 months.

Regular Set

- #7602 - Super Light
- #5030 - Light

2 - 50ml cartridges
14 - each Mixing/IntraOral tips

- #5020 - Medium
- #5010 - Heavy

2 - 50ml cartridges
10 - Mixing tips



Super Light

Shelf Life: 36 months.

Fast Set

- #7601 - Super Light
- #5031 - Light

2 - 50ml cartridges
14 - each Mixing/IntraOral tips

- #5021 - Medium
- #5011 - Heavy

2 - 50ml cartridges
10 - Mixing tips



Light

Shelf Life: 36 months.

Regular Set

- #5032 - Light

10 - 50ml cartridges
70 - each Mixing/IntraOral tips

- #5022 - Medium
- #5012 - Heavy

10 - 50ml cartridges
50 - Mixing tips



Medium

Shelf Life: 36 months.

Fast Set

- #5033 - Light

10 - 50ml cartridges
70 - each Mixing/IntraOral tips

- #5023 - Medium
- #5013 - Heavy

10 - 50ml cartridges
50 - Mixing tips



Heavy

Shelf Life: 36 months.

Clone™ Bite

BITE REGISTRATION MATERIAL

- Exact bite registration
- Fast 60-second set
- Thixotropic
- Very hard: 95 Shore A



Clone Bite bite registration material is perfect for preliminary impressions for crown/bridge procedures and temporaries.



Clone Bite makes an accurate impression for an exact bite registration.

#7600 - Clone Bite

2 - 50ml cartridge
12 - Mixing tips

Shelf Life: 36 months.



Chromaclone® PVS Accessories

#5040 - Chromaclone PVS Putty

1 - each 250g Base and 250g catalyst
2 - Scoops



#5100 - Dispensing Gun

1 - Gun



#5101 - Tray Adhesive

1 - 10ml Bottle



Chromaclone® PVS Mixing Tips

- Auto-mixing
- Color-coded to match different consistencies

Super Light = Yellow Mixing Tip

Light = Yellow Mixing Tip

Medium = Pink Mixing Tip

Heavy = Green Mixing Tip

CloneBite = Green Mixing Tip

Yellow IntraOral Tips attach to Yellow Mixing Tips

#2902 - Yellow

#2903 - Pink

#2904 - Green

50 - Tips



#2906 - IntraOral Impression

50 - Tips





Opalescence®
trèswhite®
..... **SUPREME**

TrèsFast & Dazzling

15% hydrogen peroxide
for on-the-go whitening.

new
15%!

